

Unicare Black Latex (GS290*)

black · medical grade · medium weight latex · examination glove

Physical Specifications:

Material: Latex
Colour: Black
Length: 240 mm
Weight: 5.5 g
Thickness: 0.08 mm
Finish: Powder free
Fit: Ambidextrous
Sterile: No

Order Codes and Size Information:

GS2902 | Small
GS2903 | Medium
GS2904 | Large
GS2905 | Extra large

Packaging:

Each dispenser box contains 100 gloves
Each case contains 10 dispenser boxes

Classification:



For more about the Unicare range, please click on the image



Features


- AQL 1.5 and EN455 parts 1-4
- Tested to EN374 Type B
- Free from powder
- Free from phthalates
- Tested to EN1186
- Beaded cuff
- Manufactured with textured palms and fingertips
- Medium-weight construction
- Touchscreen compatible

Benefits

- Meets medical standards including pinholes
- Resistant to a range of chemicals
- Reduces risks of latex allergies (Type I sensitisation) and minimises particulate contamination
- Ambidextrous fit means you don't have to split a pair to replace a damaged glove
- Suitable for contact with all food types
- Beaded cuff provides extra strength during donning and helps prevent liquid roll back
- Textured palms and fingertips provide improved grip, particularly when working in wet or oily conditions
- Medium-weight construction provides a balance of strength and sensitivity
- Suitable for use with mobile and tablet displays without removing gloves

Sectors

- Tattoo
- Piercing
- Hair and beauty
- Cleaning and janitorial
- Emergency services
- Food processing
- General assembly
- Healthcare
- Janitorial
- Local Authority work
- Painting and decorating
- Printing

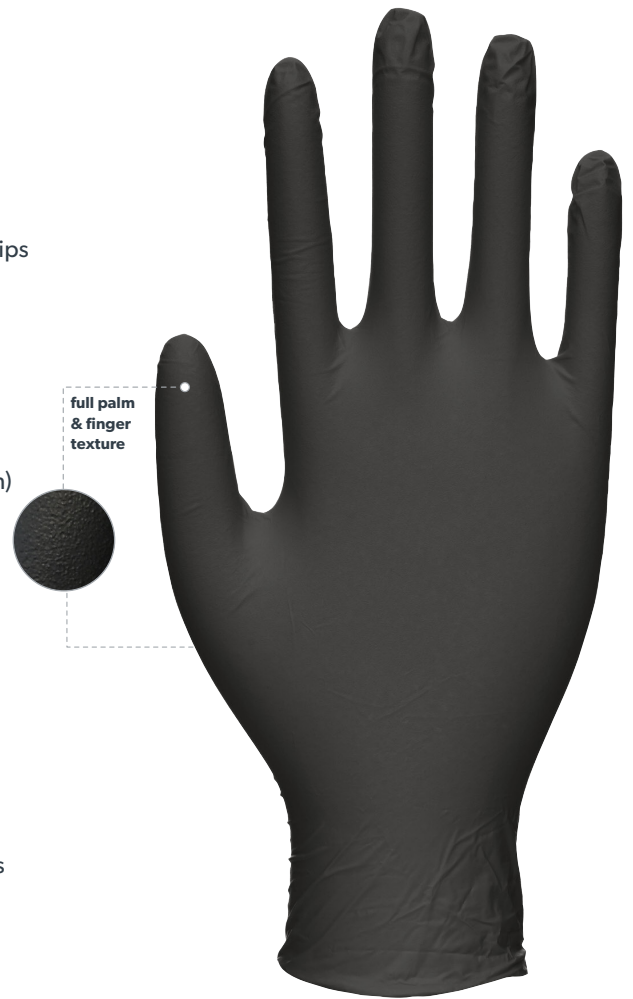

UNIGLOVES UK
Registered in England No: 4010200
VAT registration 791817787

3 Ambley Green,
Gillingham Business Park,
Gillingham, ME8 0NJ, UK
+44 (0)800 049 6602
enquiries@unigloves.co.uk
uk.unigloves.com

UNIGLOVES GERMANY
Arzt- & Klinikbedarf
Handelsgesellschaft mbH

Camp-Spich-Straße 71
D-53842 Troisdorf,
Germany
+49 (0) 2241-9323-0
unigloves@unigloves.de
unigloves.de

