

# Zero (GM005\*)

zinc, sulphur and accelerator free · low modulus · anti-static  
 blue · medical grade · medium weight · nitrile · examination gloves

## Physical Specifications:

Material: Zinc free, sulphur free, and accelerator free nitrile (latex free)  
 Colour: Blue  
 Length: 240 mm  
 Weight: 3.8 g  
 Thickness: 0.08 mm  
 Finish: Powder free, chlorinated  
 Fit: Ambidextrous  
 Sterile: No

## Order Codes and Size Information:

GM0051 | Extra small  
 GM0052 | Small  
 GM0053 | Medium  
 GM0054 | Large  
 GM0055 | Extra large

## Packaging:

Each dispenser box contains 100 gloves  
 Each case contains 10 dispenser boxes

## Classification:



For more about the **Advanced** range, please click on the image



## Features


- Free from zinc, sulphur and accelerators
- Tested to EN374 Type B
- Low modulus formulation
- Anti-static / EN1149 compliant
- AQL 1.5 and EN455 parts 1-4
- Free from powder
- Free from phthalates
- Tested to EN1186
- Nitrile formulation conforms to the hand during use
- Manufactured with textured fingertips
- Medium-weight construction
- Touchscreen compatible

## Benefits

- ZAS free formulation reduces risk of Type I and Type IV allergic reactions
- Resistant to a range of chemicals
- Low modulus formulation closely mimics the elasticity and comfort of a latex glove which minimises hand fatigue but retains the other benefits of nitrile
- Anti-static formulation suppresses static charge, preventing sparks
- Meets medical standards including pinholes
- Suitable for those with latex allergies (Type I sensitisation)
- Reduces risks of latex allergies (Type I sensitisation) and minimises particulate contamination
- Suitable for contact with all food types
- Nitrile formulation conforms to the hand during use for a more comfortable fit
- Textured fingertips provide improved grip, particularly when working in wet or oily conditions
- Medium-weight construction provides a balance of strength and sensitivity
- Suitable for use with mobile and tablet displays without removing gloves

## Sectors

- Aerospace • Chemical manufacture and use • Dentistry
- Dental hygiene • Electronics • Engineering • Food processing
- Healthcare • Laboratory work • Manufacturing and assembly
- Pharmaceutical work • Printing

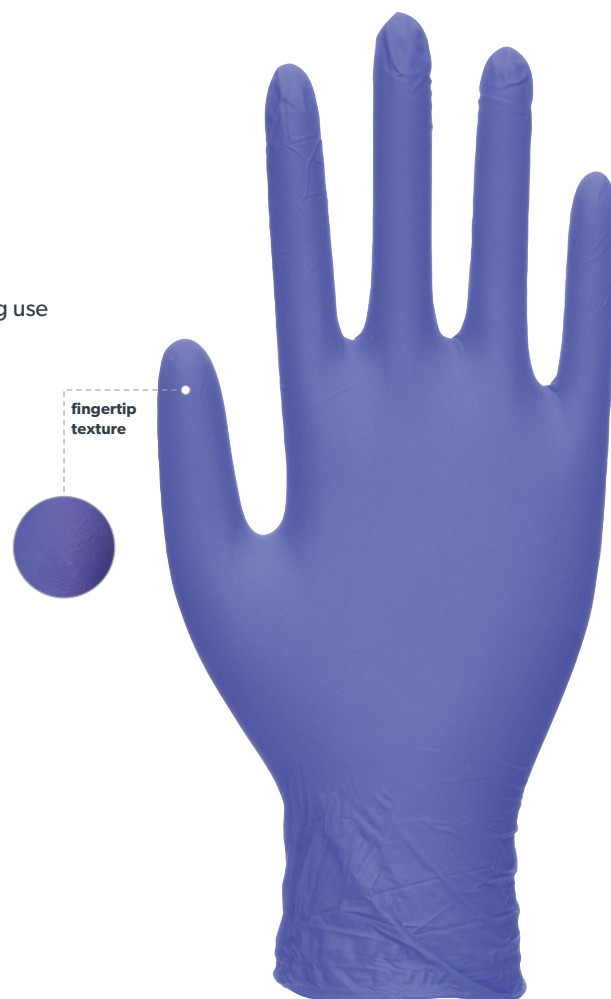

**UNIGLOVES UK**  
 Registered in England No: 4010200  
 VAT registration 791817787

3 Ambley Green,  
 Gillingham Business Park,  
 Gillingham, ME8 0NJ, UK  
 +44 (0)800 049 6602  
 enquiries@unigloves.co.uk  
 uk.unigloves.com

**UNIGLOVES GERMANY**  
 Arzt- & Klinikbedarf  
 Handelsgesellschaft mbH

Camp-Spich-Straße 71  
 D-53842 Troisdorf,  
 Germany  
 +49 (0) 2241-9323-0  
 unigloves@unigloves.de  
 de.unigloves.com

