

Unicare Safetec Latex (GS001*)

natural · medical grade · medium weight
 latex · examination gloves

Physical Specifications:

Material: Latex
 Colour: Natural
 Length: 240 mm
 Weight: 5.7 g
 Thickness: 0.11 mm
 Finish: Powder free,
 Fit: Ambidextrous
 Sterile: No

Order Codes and Size Information:

GS0011 | Extra small
 GS0012 | Small
 GS0013 | Medium
 GS0014 | Large
 GS0015 | Extra large

Packaging:

Each dispenser box
 contains 100 gloves
 Each case contains 10
 dispenser boxes

Classification:



For more about
 the Unicare
 range, please
 click
 on the
 image



Features


- AQL 1.5 and EN455 parts 1-4
- Tested to EN374 Type B
- Free from powder
- Free from phthalates
- Tested to EN1186
- Beaded cuff
- Manufactured with textured palms and fingertips
- Medium-weight construction
- Touchscreen compatible

Benefits

- Meets medical standards including pinholes
- Resistant to a range of chemicals
- Reduces risks of latex allergies (Type I sensitisation) and minimises particulate contamination
- Suitable for contact with food
- Suitable for contact with all food types
- Beaded cuff provides extra strength during donning and helps prevent liquid roll back
- Textured palms and fingertips provide improved grip, particularly when working in wet or oily conditions
- Medium-weight construction provides a balance of strength and sensitivity
- Suitable for use with mobile and tablet displays without removing gloves

Sectors

- Care home work
- Cleaning and janitorial
- Dentistry • Dental hygiene
- Emergency services
- Food processing
- General assembly
- Healthcare • Hospital and clinical duties
- Janitorial • Local Authority work
- Painting and decorating
- Printing

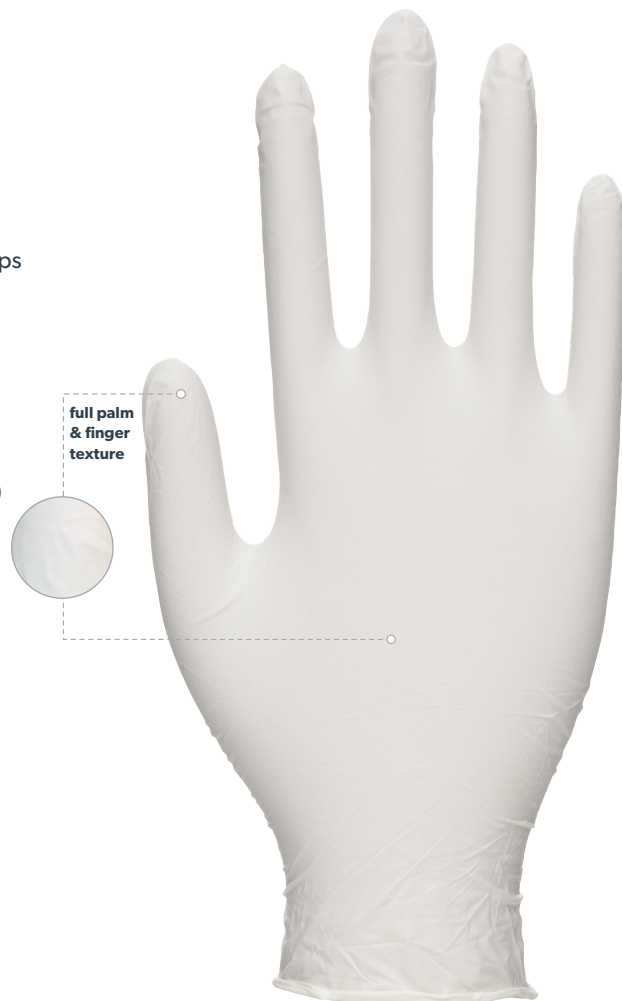

UNIGLOVES UK
 Registered in England No: 4010200
 VAT registration 791817787

3 Ambley Green,
 Gillingham Business Park,
 Gillingham, ME8 0NJ, UK
 +44 (0)800 049 6602
 enquiries@unigloves.co.uk
 uk.unigloves.com

UNIGLOVES GERMANY
 Arzt- & Klinikbedarf
 Handelsgesellschaft mbH

Camp-Spich-Straße 71
 D-53842 Troisdorf,
 Germany
 +49 (0) 2241-9323-0
 unigloves@unigloves.de
 unigloves.de



Issue No. 6

Prepared at Gillingham
 on 09/06/25