



LANO-E-GEL (190X*)

moisturising · offline washed · natural · medical grade
 medium weight · latex · examination gloves

Physical Specifications:

Material: Premium grade
 Latex (offline chlorinated)
 Colour: Natural
 Length: 240 mm
 Weight: 6.0 g
 Thickness: 0.13 mm
 Finish: Powder free,
 chlorinated
 Fit: Ambidextrous
 Sterile: No

Order Codes and Size Information:

1901 | Extra small
 1902 | Small
 1903 | Medium
 1904 | Large
 1905 | Extra large

Packaging:

Each dispenser box
 contains 100 gloves
 Each case contains 10
 dispenser boxes

Classification:



For more about
 the **Premium**
 range, please
 click
 on the
 image



Features


- Washed offline
- Microfilm inner coating contains a combination of natural moisturisers
- Tested to EN374
- AQL 1.5 and EN455 parts 1-4
- Free from powder
- Tested to EN1186
- Beaded cuff
- Manufactured with textured palms and fingertips
- Medium-weight construction
- Touchscreen compatible

Benefits

- Offline washing ensures ease of donning, increased comfort in use and reduced risk of allergic reactions
- Microfilm inner coating helps protect and nourish skin, and aids donning
- Resistant to a range of chemicals
- Meets medical standards including pinholes
- Reduces risks of latex allergies (Type I sensitisation) and minimises particulate contamination
- Suitable for contact with all food types
- Beaded cuff provides extra strength during donning and helps prevent liquid roll back
- Textured palms and fingertips provide improved grip, particularly when working in wet or oily conditions
- Medium-weight construction provides a balance of strength and sensitivity
- Suitable for use with mobile and tablet displays without removing gloves

Sectors

- Dentistry • Dental hygiene
- Floristry and gardening
- Hair and beauty
- Healthcare

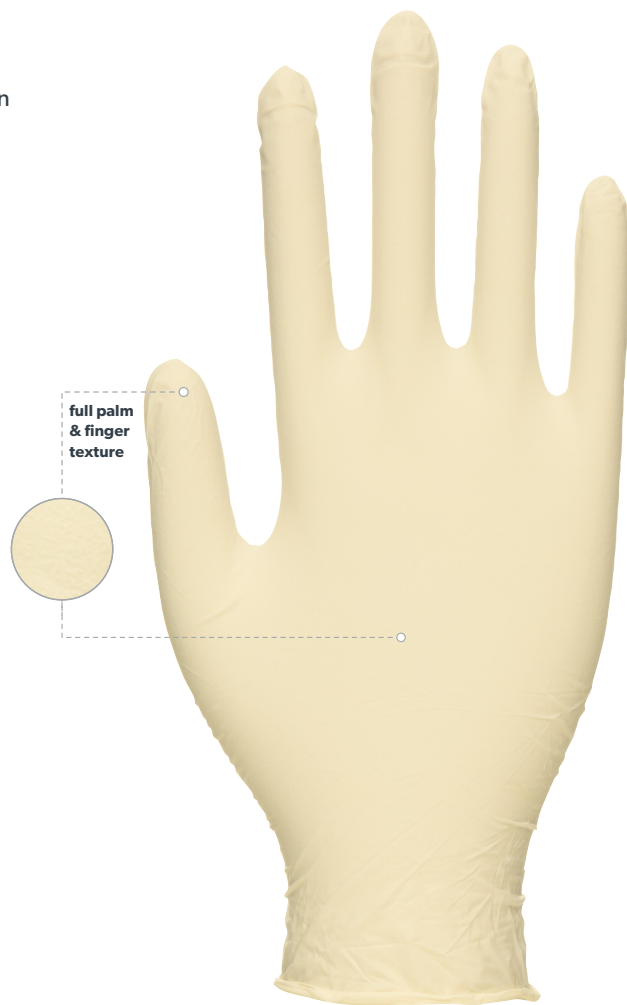

UNIGLOVES UK
 Registered in England No: 4010200
 VAT registration 791817787

3 Ambley Green,
 Gillingham Business Park,
 Gillingham, ME8 0NJ, UK
 +44 (0)800 049 6602
 enquiries@unigloves.co.uk
 uk.unigloves.com

UNIGLOVES GERMANY
 Arzt- & Klinikbedarf
 Handelsgesellschaft mbH

Camp-Spich-Straße 71
 D-53842 Troisdorf,
 Germany
 +49 (0) 2241-9323-0
 unigloves@unigloves.de
 de.unigloves.com



Issue No. 1

Prepared at Gillingham
 on 29/09/25